



AFRL-HE-BR-TR-2006-0075

**ADVANCES IN BIOTECHNOLOGY AND THE BIOSCIENCES FOR
WARFIGHTER PERFORMANCE AND PROTECTION:
ANTI-APTAMERS FOR ENVENOM**

**George W. Irving
Amanda J. Tijerina
Mark Sloan**

**Conceptual MindWorks, Inc.
9830 Colonnade Blvd. Ste 377
San Antonio TX 78230**

October 2006

Approved for public release, distribution unlimited

**Air Force Research Laboratory
Human Effectiveness Directorate
Directed Energy Bioeffects Division
Counter-Proliferation Branch
Brooks City-Base TX 78235**

NOTICE

Using Government drawings, specifications, or other data included in this document for any purpose other than Government procurement does not in any way obligate the U.S. Government. The fact that the Government formulated or supplied the drawings, specifications, or other data does not license the holder or any other person or corporation; or convey any rights or permission to manufacture, use, or sell any patented invention that may relate to them.

This report was cleared for public release by the Air Force Research Laboratory, Brooks City-Base Public Affairs Office (PAO) and is releasable to the National Technical Information Service (NTIS). It will be available to the general public, including foreign nationals.

National Technical Information Service
5285 Port Royal Road, Springfield VA 22161

Federal Government agencies and their contractors registered with Defense Technical Information Center should direct requests for copies of this report to:

Defense Technical Information Center
8725 John J. Kingman Rd., STE 0944, Ft Belvoir VA 22060-6218

TECHNICAL REVIEW AND APPROVAL

AFRL-HE-BR-TR-2006-0075

**THIS TECHNICAL REPORT HAS BEEN REVIEWED AND IS APPROVED FOR
PUBLICATION.**

//SIGNED//
SHARON K. GARCIA
Contract Monitor

//SIGNED//
F. WESLEY BAUMGARDNER
Deputy Chief, Biosciences and Protection Division

This report is published in the interest of scientific and technical information exchange and its publication does not constitute the Government's approval or disapproval of its ideas or findings.

REPORT DOCUMENTATION PAGE				Form Approved OMB No. 0704-0188	
Public reporting burden for this collection of information is estimated to average 1 hour per response, including the time for reviewing instructions, searching existing data sources, gathering and maintaining the data needed, and completing and reviewing this collection of information. Send comments regarding this burden estimate or any other aspect of this collection of information, including suggestions for reducing this burden to Department of Defense, Washington Headquarters Services, Directorate for Information Operations and Reports (0704-0188), 1215 Jefferson Davis Highway, Suite 1204, Arlington, VA 22202-4302. Respondents should be aware that notwithstanding any other provision of law, no person shall be subject to any penalty for failing to comply with a collection of information if it does not display a currently valid OMB control number. PLEASE DO NOT RETURN YOUR FORM TO THE ABOVE ADDRESS.					
1. REPORT DATE (DD-MM-YYYY) 24 May 2006		2. REPORT TYPE Technical Report		April 2005 – May 2006	
4. TITLE AND SUBTITLE Advances in Biotechnology and the Biosciences for Warfighter Performance and Protection: Anti-Aptamers for Envenom				5a. CONTRACT NUMBER FA8650-05-C-6520	
				5b. GRANT NUMBER	
				5c. PROGRAM ELEMENT NUMBER 62202F	
6. AUTHOR(S) George W. Irving Amanda J. Tijerina Mark Sloan				5d. PROJECT NUMBER 7757	
				5e. TASK NUMBER P4	
				5f. WORK UNIT NUMBER 31	
7. PERFORMING ORGANIZATION NAME(S) AND ADDRESS(ES) Conceptual MindWorks, Inc. 9830 Colonnade Blvd., Suite 377 San Antonio, TX 78230				8. PERFORMING ORGANIZATION REPORT NUMBER	
Air Force Research Laboratory Human Effectiveness Directorate Biosciences and Protection Division 2486 Gillingham Dr Brooks City-Base, TX 78235-5105				10. SPONSOR/MONITOR'S ACRONYM(S) AFRL; AFRL/HEPC	
				11. SPONSOR/MONITOR'S REPORT NUMBER(S) AFRL-HE-BR-TR-2006-0075	
12. DISTRIBUTION / AVAILABILITY STATEMENT Distribution unlimited					
13. SUPPLEMENTARY NOTES Air Force Research Laboratory Program Manager: Dr Sharon K. Garcia, Technical Program Manager: Dr Johnathan Kiel					
14. ABSTRACT This effort was focused on developing a novel, aptamer-based antivenin for treatment of envenomation by the Kurdistan Viper (<i>Vipera raddei kurdistanica</i>). The research was conducted to provide evidence to prove whether a synthetic, aptamer-based antivenin could be developed to treat snake envenomations in humans. Using PLA ₂ from <i>Crotalus durissus terrificus</i> venom as a simulant of the Kurdistan viper venom, two tissue culture cell lines were examined and developed for <i>in vitro</i> cell culture models. For each cell line, an LD ₅₀ value was determined post PLA ₂ exposure at various concentrations. Cytotoxicity activity was determined by utilizing an XTT colorimetric assay. DNA aptamers developed against the PLA ₂ were tested in these <i>in vitro</i> models, along with known PLA ₂ inhibitors. Inhibitors were tested for their effectiveness against these LD ₅₀ values for each cell line. However, in these assays, known LD ₅₀ values for PLA ₂ did not prove to be toxic to the cells themselves. Higher concentrations of PLA ₂ , were also ineffective in killing cells the effectiveness of the inhibitor's ability to decrease PLA ₂ activity, thereby preventing toxicity to cells, could not be determined. Additionally, no determination was able to be made on the efficacy of the aptamers.					
15. SUBJECT TERMS Kurdistan Viper, aptamer-based antivenin; envenomation treatments					
16. SECURITY CLASSIFICATION OF:			17. LIMITATION OF ABSTRACT	18. NUMBER OF PAGES	19a. NAME OF RESPONSIBLE PERSON
a. REPORT Uncl	b. ABSTRACT Uncl	c. THIS PAGE Uncl			Sharon Garcia
			Uncl	42	19b. TELEPHONE NUMBER (include area code)

THIS PAGE INTENTIONALLY LEFT BLANK

Table of Contents

Table of Contents	i
List of Tables	iii
List of Figures	iii
List of Acronyms	iv
PREFACE	v
EXECUTIVE SUMMARY	1
CHAPTER 1 PROGRAM DESCRIPTION	
1.1 Technical Approach	2
1.2 Customer	2
1.3 Goals and Objectives of the Research and Development	3
1.4 Customer Benefits/Payoffs	3
CHAPTER 2 TECHNICAL DISCUSSION	
2.1 Background	3
2.2 Hypotheses	5
CHAPTER 3 MATERIALS AND METHODS	
3.1 Materials	6
3.1.1 Phospholipase A ₂ and Inhibitors	6
3.1.2 Cell Lines	6
3.2 Methods	8
3.2.1 Cell Subculturing	8
3.2.2 Cell Harvesting	8
3.2.3 Determining Cell Concentration	9
3.2.4 Confluency Testing	9
3.2.5 Cytotoxicity Assay	10

CHAPTER 4 RESULTS

4.1	Confluency Testing	11
4.2	LD ₅₀ Values	12
4.3	Inhibition Evaluation.....	13

CHAPTER 5 MSDS

5.1	Phospholipase A ₂ from <i>Crotalus durissus terrificus</i>	15
5.2	Diethylenetriaminepentaacetic acid	19
5.3	D,L-Erythro-dihydrosphingosine	24
5.4	7,7-Dimethyl-(5Z, 8Z)-Eicosadienoic acid	27

CHAPTER 6 REFERENCES

List of References	33
--------------------------	----

List of Tables

Table 1. Results of confluency percentages from various starting concentrations after 24, 48 and 72 hour culturing periods	12
--	----

List of Figures

Figure 1 The Kurdistan Viper (<i>Vipera raddei kurdistanica</i>)	2
Figure 2 Extensive damage to forearm after envenomation	4
Figure 3 Phase contrast micrograph of C2C12 mouse myoblast cultured cells (ATCC: CRL-1772) at low density and high density confluency	7
Figure 4 Phase contrast micrograph of MDCK cultured cells (ATCC: CCL-34) at low density and high density confluency	8
Figure 5 Experimental design fro confluency tests of C2C12 and MDCK cell lines	9
Figure 6 Experimental design of XTT assays with varying concentrations of PLA ₂	10
Figure 7 Experimental design of XTT assays with varying inhibitors of PLA ₂	11
Figure 8 Cell viability of C2C12 cells 2 and 4 hours post exposure to various concentrations of PLA ₂	12
Figure 9 Cell viability of MDCK cells 2 and 4 hours post exposure to various concentrations of PLA ₂	13
Figure 10 C2C12 viability 2 and 4 hours post exposure to PLA ₂ and it known inhibitors	14
Figure 11 MDCK viability 2 and 4 hours post exposure to PLA ₂ and its known inhibitors	14

List of Acronyms

AFSOC	Air Force Special Operations Command
AMC	Air Mobility Command
ATCC	American Type Culture Collection
CMI	Conceptual MindWorks, Inc.
DESA	7,7-Dimethyl-(5Z,8Z)-eicosaienoic acid
DHS	DL- <i>erythro</i> -Dihydrosphingosine
DMEM	Dulbecco's Modified Eagle's Medium
DTPA	Diethylenetriaminepentaacetic acid
EMEM	Eagle's Minimum Essential Medium
FBS	Fetal Bovine Serum
PLA ₂	Phospholipase A ₂
SOCOM	Special Operations Command

PREFACE

The scope of this effort was to address the development of a novel, aptamer-based anti-venom for treatment of envenomation by the Kurdistan Viper (*Vipera raddei kurdistanica*) and to provide evidence for whether or not a synthetic, aptamer-based antivenin can be developed which could be used to treat snake envenomations in humans. The development of this antivenin would be groundbreaking in the neutralization of toxins and would provide an unfilled ability to treat warfighters envenomated by this organism.

No classified information was used to generate, or is included in, this report. In accordance with AFRL program guidance, this report was submitted to the project sponsor for further disposition. Any restrictions upon the dissemination of this report, or any classification of the information contained herein, are at the discretion of the United States Air Force, Counterproliferation Branch, Brooks City-Base, Texas.

EXECUTIVE SUMMARY

This effort was focused on developing a novel, aptamer-based antivenin for treatment of envenomation by the Kurdistan Viper (*Vipera raddei kurdistanica*). The research was conducted to provide evidence to prove whether a synthetic, aptamer-based antivenin could be developed to treat snake envenomations in humans. The development of this antivenin would be groundbreaking in the neutralization of toxins and therefore provide an unfilled ability to treat warfighters envenomated by this organism.

During this 1 year effort, the Conceptual MindWorks (CMI) team was able to address two of the specific aims of the project. Using PLA₂ from *Crotalus durissus terrificus* venom as a simulant of the Kurdistan viper venom (no source available at this time), two tissue culture cell lines were examined and developed for *in vitro* cell culture models. The two cell lines used for this study were C2C12 (CRL-1772), a mouse myoblast cell line, and MDCK (CCL-34), a canine kidney cell line. For each cell line, an LD₅₀ value was determined post PLA₂ exposure at various concentrations. Cytotoxicity activity was determined by utilizing an XTT colorimetric assay.

DNA aptamers developed against the PLA₂ were tested in these *in vitro* models, along with known PLA₂ inhibitors. Inhibitors were tested for their effectiveness against these LD₅₀ values for each cell line. However, in these assays, known LD₅₀ values for PLA₂ did not prove to be toxic to the cells themselves. Higher concentrations of PLA₂, 10µg/mL for C2C12 cells and 50µg/mL for MDCK cells, were also ineffective in killing cells the effectiveness of the inhibitor's ability to decrease PLA₂ activity, thereby preventing toxicity to cells, could not be determined. Additionally, no determination was able to be made on the efficacy of the aptamers.

This research was designed to extend the research, planning and integration efforts underway by the Air Force Research Laboratory (AFRL) to produce techniques for the neutralization of biowarfare agents.

Chapter 1. Program Description

1.1 Technical Approach

A major component found in the venom of the Kurdistan Viper, shown in Figure 1, known to cause altered pathophysiological effects associated with post-envenomation, such as hemorrhage and necrosis, is Phospholipase A₂ (PLA₂).^{1,2}



Figure 1. The Kurdistan Viper (*Vipera raddei kurdistanica*)

CMI's research team will develop highly specific and selective aptamers against PLA₂, found in the venom from the viper, using a SELEX like process. The team will determine which mouse and/or human cell lines are susceptible. Possible cell lines are from mouse muscle tissue as the final anti-PLA₂ testing will be in a mouse model. Using an *in vitro* culture model, the team will select aptamers with the greatest ability to neutralize PLA₂. Cytotoxicity will be measured using two methods: 1) enzymatic assays and 2) neutral red viability test. This will allow quantification of cell viability. Variables will be dose of PLA₂ and concentration of aptamer cocktail. During the third year of this research, the aptamer production and cell culture work will be extended to the mouse model. BALB/c mice will be used for *in vivo* assays; assuming that the *in vitro* assays demonstrate that an anti-PLA₂ toxins aptamer cocktail spares cells from the applied venom. Protection assays will involve post-treatment of mice with aptamer cocktail twice-daily by subcutaneous injections.³

1.2 Customer

The primary customer is the Air Force Special Operations Command (AFSOC) and US Special Operations Command (USSOCOM). In general, the Departments of Defense, Homeland Security, and Environmental Protection Agency will also benefit from this capability.

1.3 Goals and Objectives of the Research and Development

Specific aims of the project are as follows:

1. Development of DNA based aptamers against PLA₂ found in the venom from Kurdistan Viper.
2. Selection of DNA aptamers that provide the greatest protection from the pathogenic effects of PLA₂ *in vitro* cell culture model.
3. Development of a mouse model for poisonous snake envenomation.
4. Testing of aptamer cocktail in a mouse model to determine if DNA based aptamers given post exposure to PLA₂ spares the mice.

1.4 Customer Benefits/Payoffs

Testing the efficacy of aptamers selected against PLA₂ will be carried out using appropriate cell lines and a range of cytotoxicity, cell proliferation and function tests. Additional benefits of this research will be knowledge gained for similar studies on excreted toxins from bacterial select agents. A further benefit of considering the use of derivative dendrite polymers or aptamers in any studies is the possibility that these could be designed to bind to a wide range of bacteria and viruses and also block their activity.⁴

Chapter 2. Technical Discussion

2.1 Background

Venomous snakes represent a potentially serious threat to US and coalition forces deployed to Iraq. The risk of venomous snakebites will be higher for ground forces during field operations. Although total numbers of bite casualties will likely be very small, victims are at risk for sustaining serious injury, disability and possible death. Successfully treating these casualties requires knowledge of emergency first aid and management of the bite site. The most effective treatment for significant snake envenomation (the act of injecting venom by a bite) is geographic and species specific antivenin administered by medical personnel. The most important decision will likely be determining if antivenin should be administered; not all snakebites result in significant envenomation requiring antivenin. The degree of envenomation is judged according to clinical criteria such as the presence of widely distributed pain, edema progressing toward the trunk, petechiae (pinpoint reddish rash) or ecchymosis (hemorrhagic spots), and systemic symptoms including fever, nausea or vomiting. Clinicians should be aware that most antivenins are produced from horses and approximately 15

to 20 percent of patients receiving equine-based antivenins will exhibit adverse side effects⁵.

Most antivenin is produced by injecting snake venom in increasingly higher doses into horses, thereby inducing the animal's immune system to produce antibodies to the venom. The horse's blood is collected and processed to manufacture antivenin, a process that can take 6 to 8 weeks. Antivenin is marketed freeze-dried in glass vials or as a liquid in ampoules. Most antivenin had a shelf-life of approximately 3 years.⁶

The quality of antivenin and efficiency of production are directly related to the age and health of both the biting snakes and horses used in manufacturing the antivenin as well as quality control practices in the facility that produces the antivenin. The World Health Organization has published guidelines for properly producing antivenin and recommends using horses between ages 5 and 10 to produce antivenin. Approximately 15 to 20 percent of recipients of horse-derived antivenin may demonstrate early adverse pyrogenic, anaphylactic or anaphylactoid reactions to current snake antivenins. Medical personnel administering antivenin should be trained and equipped to treat adverse antivenin reactions.⁷

Military personnel could sustain a wound from a snake for which no specific antivenin is available. In these instances, one option is using polyvalent antivenin that might contain genus-specific antivenin. Six of Iraq's venomous snakes are classified as true vipers. These snakes produce hemotoxic venom that causes severe damage to blood cells and tissue of



Figure 2. Extensive damage to forearm after envenomation.

bite victim,⁸ as shown in Figure 2.

In recent years several techniques have been devised for the synthesis and functional screening of large numbers of organic molecules. Libraries of peptides, antibodies or partially randomized proteins displayed on the surfaces of cells and variants with desirable binding or catalytic abilities can be iteratively selected and amplified.^{9,10} Aptamers are artificial nucleic acid

ligands that can be generated against amino acids, drugs, proteins and other molecules. They are isolated from complex libraries of synthetic nucleic acids by an iterative process of adsorption, recovery and reamplification, a process known as SELEX (Systematic Evolution of Ligands by Exponential Enrichment).^{11,12} Today, the SELEX process has been applied to more than a hundred different target molecules. A wide variety of molecules have been targeted by *in vitro* selection experiments that have yielded specific nucleic acid-binding sequences. Perhaps the smallest target being chelated zinc molecule (atom).¹³ Literally hundreds of proteins have been targeted by SELEX selection and yielded anti-protein aptamers.^{14,15} These proteins include toxins,^{16,17} glycoproteins,¹⁸ HIV 1 Rev,¹⁹ immunoglobulins, signal molecules,²⁰ growth factors,²¹ and toxins. Aptamers have been produced against very diverse molecules ranging from organic dyes to nucleotides, to amino acids, peptides and complex proteins.²²

Furthermore, at Brooks City-Base, our own laboratory selected aptamers against live anthrax spores,^{17,23} shiga toxin, whole virus VEE and whole bacterial cells showing that aptamers can be obtained for almost any desired target whether simple or complex.

2.2 Hypotheses

Antivenin therapy is the mainstay of medical treatment of snakebites, along with administration of plasma expanders, pain medication, diazepam, tetanus toxoid, antiseptics and antibiotics. Patients who have pain, swelling, ecchymosis, systemic symptoms or abnormal laboratory findings within 30 minutes to one hour of a bite are probable candidates to receive antivenin therapy. Before receiving antivenin therapy, the patient must be tested for hypersensitivity to the antivenin. Antivenin therapy is the most effective when given within four hours of a snakebite.²⁴

Hemotoxic venom attacks the circulatory system and muscle tissue causing excessive scarring, gangrene, permanent disuse of motor skills and sometimes leads to amputation of the affected area.

An anti-PLA₂ aptamer could fulfill almost all the criteria that describe an ideal antivenin: economically affordable; possess a long shelf-life under various storage conditions; exact affinity for the specific venom selected against; zero allergenicity; synthetically manufactured.

The high affinity and specificity of nucleic acid aptamers, their lack of immunogenicity and the ability to raise aptamers against any target for which an *in vitro* selection method can be devised, makes them tempting candidates for drug discovery. In the area of potential therapeutic aptamer antagonists of the toxin ricin have been isolated with IC₅₀ values in the nanomolar range.¹⁶ Furthermore, a PEG-conjugated aptamer to vascular endothelial growth factor that inhibits pathogenic angiogenesis has already

reached clinical trials and is intended for the treatment of blindness induced by macular degeneration. It was also found to be safe and relatively long lasting following injection and has a half-life in plasma of little over 9 hours^{25, 26}.

With this background we will test the following hypotheses:

1. Anti-Kurdistan Viper PLA₂ (aKVPLA₂) aptamer cocktails would provide cyto-protection to Mouse C2C12 muscle cells and MDCK canine kidney cells *in vitro* against envenomation from Kurdistan viper venom toxicity.
2. An aptamer cocktail could act as a synthetic antivenin, protecting mice from snake envenomation toxicity.

Chapter 3 Materials and Methods

3.1 Materials

3.1.1 Phospholipase A₂ and Inhibitors

PLA₂ from *Crotalus durissus terrificus* venom (P5910) and its known inhibitors: Diethylenetriaminepentaacetic acid (DTPA D6518), DL-erythro-Dihydrosphingosine (DHS D6908) and 7, 7-Dimethyl-(5Z, 8Z)-eicosaienoic acid (DESA D8008) were purchased from Sigma-Aldrich. PLA₂, DHS and DESA were reconstituted in dH₂O to make a 5mg/mL, 83mM and 30mM concentration, respectively. The stock concentrations were stored in 50μL aliquots at -20°C. DTPA was stored at room temperature and reconstituted to make appropriate stock concentrations for each individual experiment.

3.1.2 Cell Lines

The two cell lines used for this study are shown in Figures 3 and 4. C2C12 (CRL-1772), a mouse myoblast cell line, and MDCK (CCL-34), a canine kidney cell line, were both purchased from American Type Culture Collection (ATCC). Cell lines C2C12 and MDCK, were grown in Dulbecco's Modified Eagle's Medium (DMEM, ATCC 30-2002) and Eagle's Minimum Essential Medium (EMEM, ATCC 30-2003), respectively. Both media were supplemented with 10% Fetal Bovine Serum (FBS, Atlanta Biologicals S11150) and 5% Penicillin-Streptomycin-Glutamine 100X (Pen/Strep/Glu, Invitrogen 10378-016). In all protocols, respective media for each cell line was used and cells were incubated at 37°C with an air atmosphere of 5% CO₂. Cells were grown in 150cm² tissue culture flasks (Corning 430825). At near confluency, cells were either subcultured or harvested.

ATCC Number: **CRL-1772**
Designation: **C2C12**

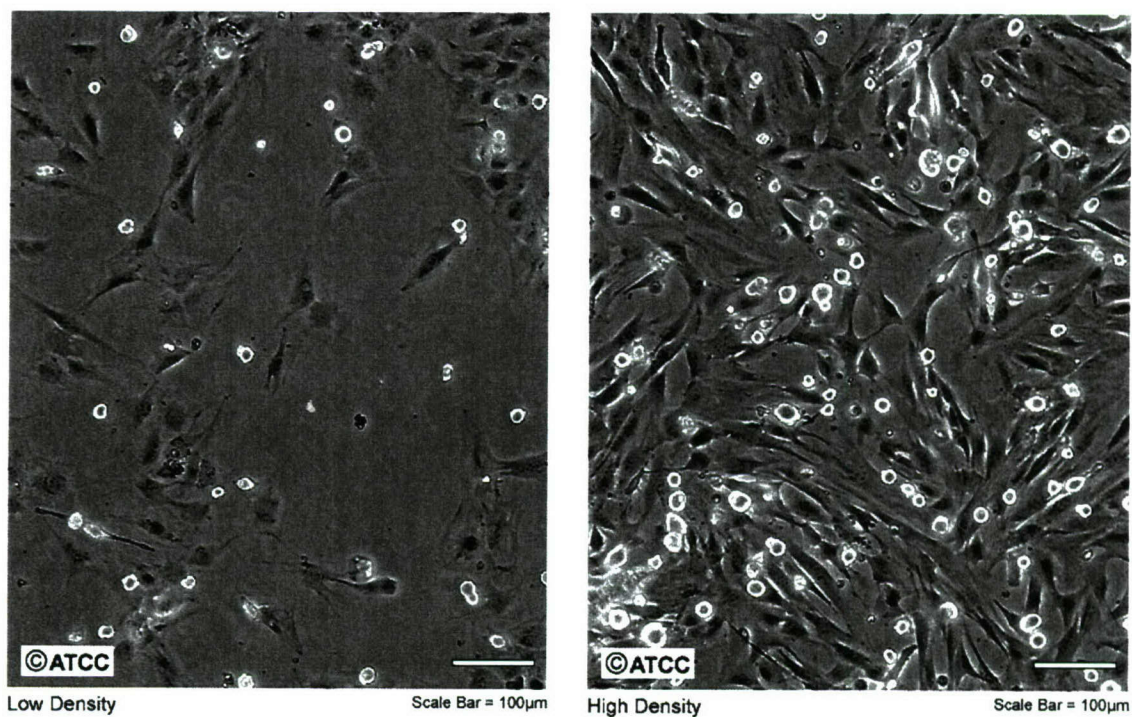


Figure 3. Phase contrast micrograph of C2C12 mouse myoblast cultured cells (ATCC: CRL-1772) at low density and high density confluency.

ATCC Number: **CCL-34**
Designation: **MDCK (NBL 2)**

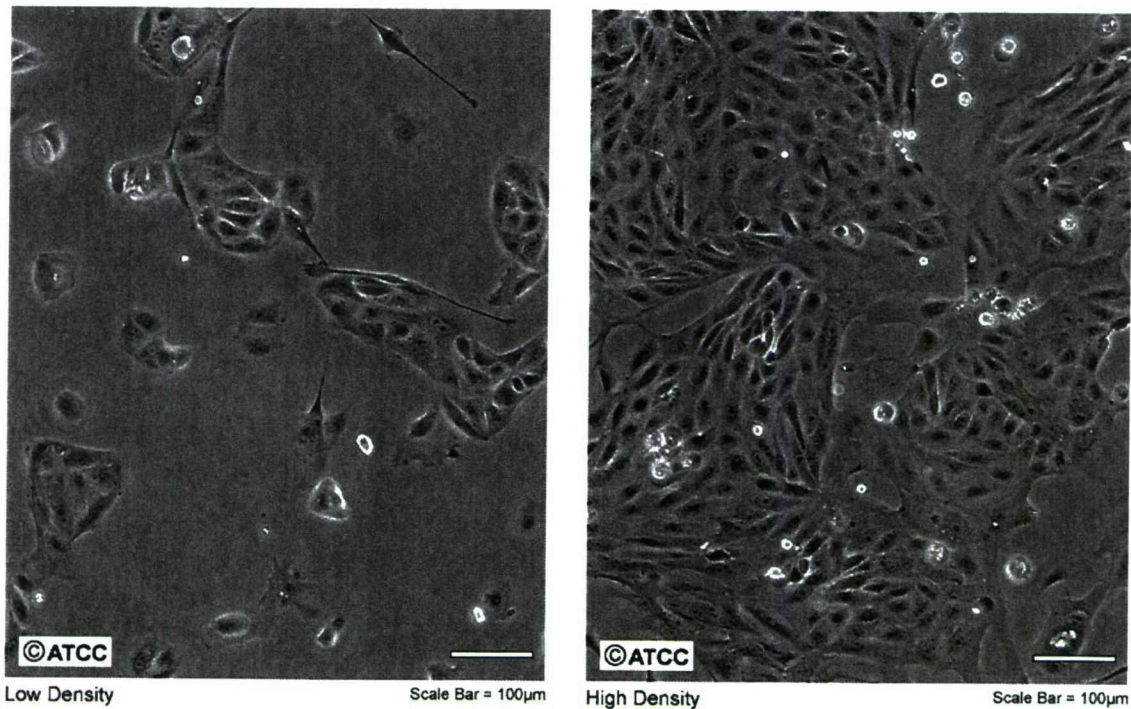


Figure 4. Phase contrast micrograph of MDCK cultured cells (ATCC: CCL-34) at low density and high density confluency.

3.2 Methods

3.2.1 Cell Subculturing

To subculture C2C12 cells, media was replaced with 6mL of Trypsin-EDTA Solution 1X (ATCC 30-2101) to resuspend cells. Cells were then incubated for 5 min. Due to the tendency of multi-layer growth in MDCK cells, 10mL of Trypsin-EDTA was used to resuspend cells for 10-15 min. Once cells were observed to have been resuspended under an inverted microscope (Nikon, TMS-F), media was added to inhibit the trypsin and give a 1:4 to 1:10 subcultivation ratio for C2C12 cells and a 1:2 to 1:6 subcultivation ratio for MDCK cells.

3.2.2 Cell Harvesting

This protocol parallels that of *cell subculturing* with the following modifications: Once cells were observed to have been resuspended under an inverted microscope, an equivalent amount of media was added to inhibit the trypsin activity. After rinsing the cells off the flask with the added media, the mixture from all respective flasks was collected in a 50mL conical tube

(Fisher 14-959-49A). The tube was centrifuged at 15°C, 3000g for 5 min. The supernatant was properly discarded and 1mL of media was added to the cells to determine cell concentration.

3.2.3 ***Determining Cell Concentration***

To determine cell concentration of a cell line, 900µL of media and 100µL of the 1mL cell/media mixture, previously mentioned, was added to a 15mL conical tube (Fisher 14-959-70C). 10µL was then pipetted into a hemocytometer (Fisher 0267110). Each grid of the hemocytometer represents a total volume of 0.1 mm³ or 10⁻⁴ cm³. Since 1 cm³ is equivalent to approximately 1 mL, the total number of cells per mL was determined using the following calculations: Cells/mL = average cell count per grid x dilution factor x 10⁴. This value was used to determine starting concentrations for confluency testing and cytotoxicity assays.

3.2.4 ***Confluency Testing***

Confluency tests were conducted for both cell lines on separate 24 well Falcon plates (Fisher 353047). Using the value of cell concentration, previously mentioned above, starting concentrations of 5.0³, 1.0⁴, 2.0⁴, and 5.0⁴ cells/mL were placed into six designated wells of each plate as shown in Figure 5 and incubated at 37°C with an air atmosphere of 5% CO₂.

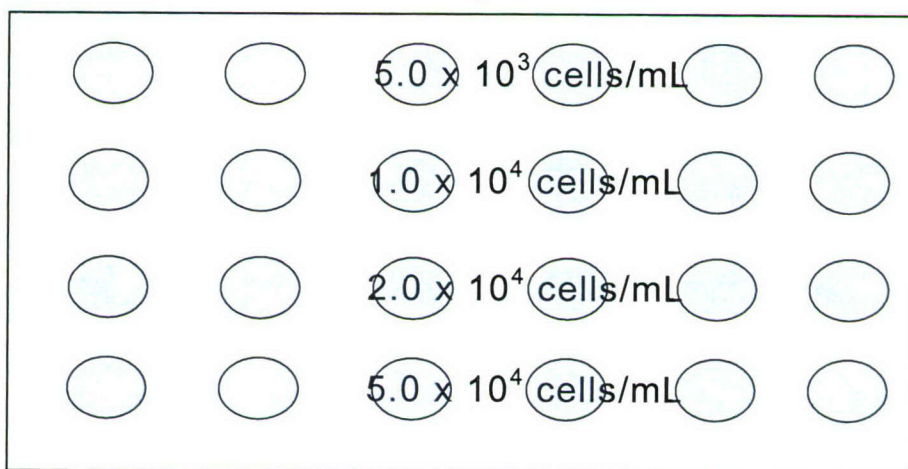


Figure 5. Experimental design for confluency tests of C2C12 and MDCK cell lines.

After incubation periods of 24- and 48-hours, each well was observed for confluency percentage on an inverted microscope. In the case of MDCK cells, an additional incubation period of 72-hours was observed. The average percentage of respective concentrations for each time point was used in determining starting concentrations for the cytotoxicity assays which followed.

3.2.5 Cytotoxicity Assay

Cytotoxicity activity was determined by utilizing an XTT Cell Proliferation Kit II (Roche 11 465 015 001). In brief, this assay is based on the cleavage of the yellow tetrazolium salt XTT to form an orange formazan dye by metabolic active cells.²⁷ Therefore, this conversion only occurs in viable cells.

When cells were at near confluency, 200 μ L of varying PLA₂ concentrations, ranging from 3 to 60 μ L/mL, were pipetted into designated wells as shown in Figure 6. Each concentration was calculated from the 5mg/mL stock and diluted with media. This assay required incubation periods of 30 min, 2, 4, and 8 hours at 37°C with an air atmosphere of 5% CO₂. One 24-well plate was used for each incubation period for each cell line.

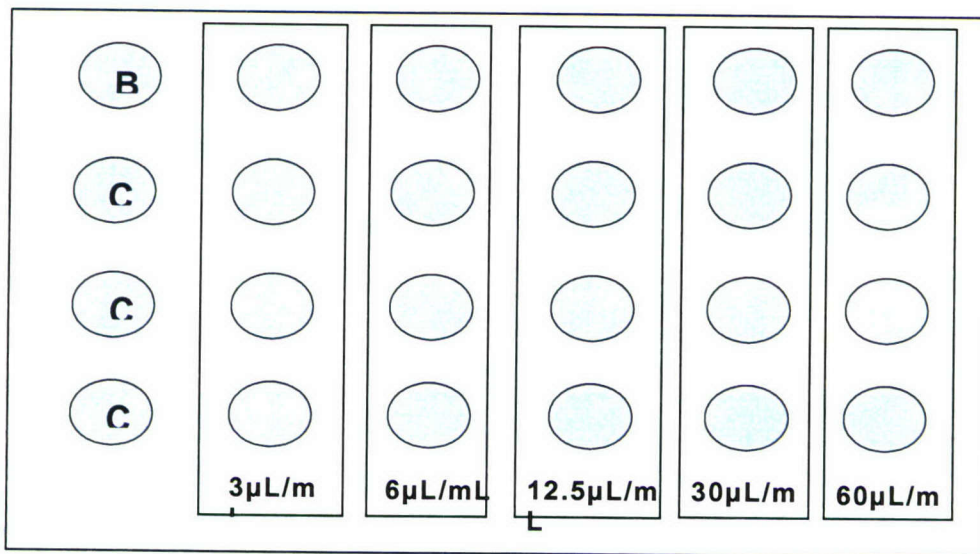


Figure 6. Experimental design of XTT assays with varying concentrations of PLA₂. (B) blank well containing only media, (C) control wells containing cells and media.

Varying concentrations of each PLA₂ inhibitor were also exposed to the cell lines. Each concentration was calculated from the 5mg/mL stock and diluted with media. Due to results demonstrating the LD₅₀ for PLA₂ when exposed to cells, this assay required incubation periods of 2 and 4 hours at 37°C with an air atmosphere of 5% CO₂. One 24-well plate was used for each incubation period for each cell line. Figure 7 shows the experimental design. When combining PLA₂ and inhibitor, 100 μ L of each concentration for

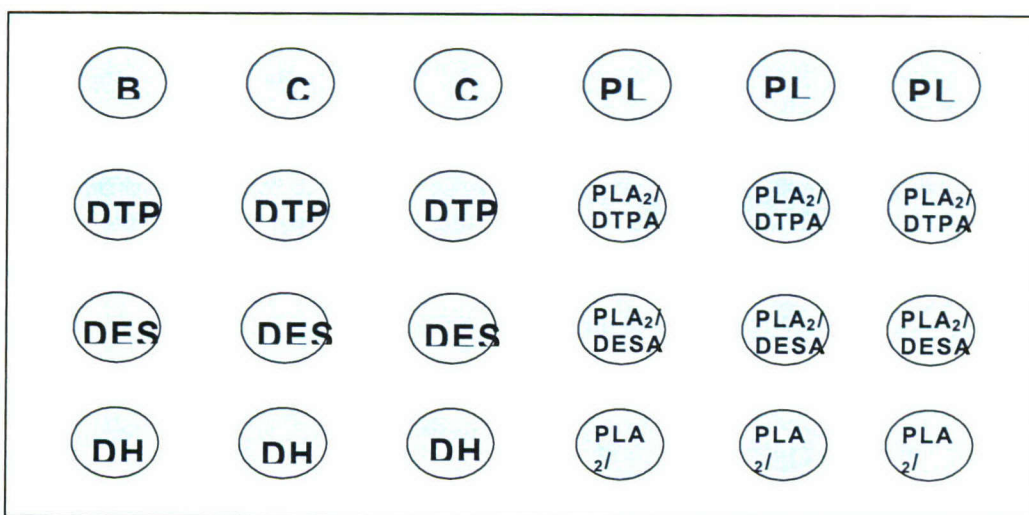


Figure 7. Experimental design of XTT assays with varying inhibitors of PLA₂. blank (B) well containing only media, control (C) wells containing cells and media, PLA₂ (PLA₂), Diethylenetriaminepentaacetic acid (DTPA), 7, 7-Dimethyl-(5Z, 8Z)-eicosaienoic acid (DESA) and DL-erythro-Dihydrosphingosine (DHS).

After incubation periods, 100µL cocktail mix of 2.8mL labeling reagent and 56µL electron-coupling reagent was pipetted into each well. The plate was then incubated for 30 minutes prior to being read by a microplate reader set at 490nm (BioTek, Synergy™ HT).

Chapter 4 Results

4.1 Confluency Testing

Table 1 shows the confluency rate of both C2C12 and MDCK cell lines. As shown, MDCK showed to have a lower proliferation rate than that of C2C12 cells. These results helped determine a set time frame for which the researcher could conduct XTT assays.

MDCK CONFLUENCY			Starting Concentration	C2C12 CONFLUENCY	
24 Hr	48 Hr	72 Hr		24 Hr	48 Hr
10%	15%	20%	5.0×10^3	30%	45%
20%	25%	30%	1.0×10^4	50%	70%
40%	70%	75%	2.0×10^4	70%	90%
50%	80%	90%	5.0×10^4	90%	100%

Table 1. Results of confluency percentages from various starting concentrations after 24, 48 and 72 hour culturing periods.

4.2 LD₅₀ Values

For each cell line, an LD₅₀ value was determined post PLA₂ exposure at various concentrations. With regard to C2C12 cells, as seen in Figure 8, the LD₅₀ value is approximately 5µg/mL for 2 hour exposure periods and approximately 4µg/mL for 4 hour exposure periods. Interestingly, an increase in cell proliferation between 30 minute and 2 hours exposure periods occurred regularly in C2C12 XTT assays.

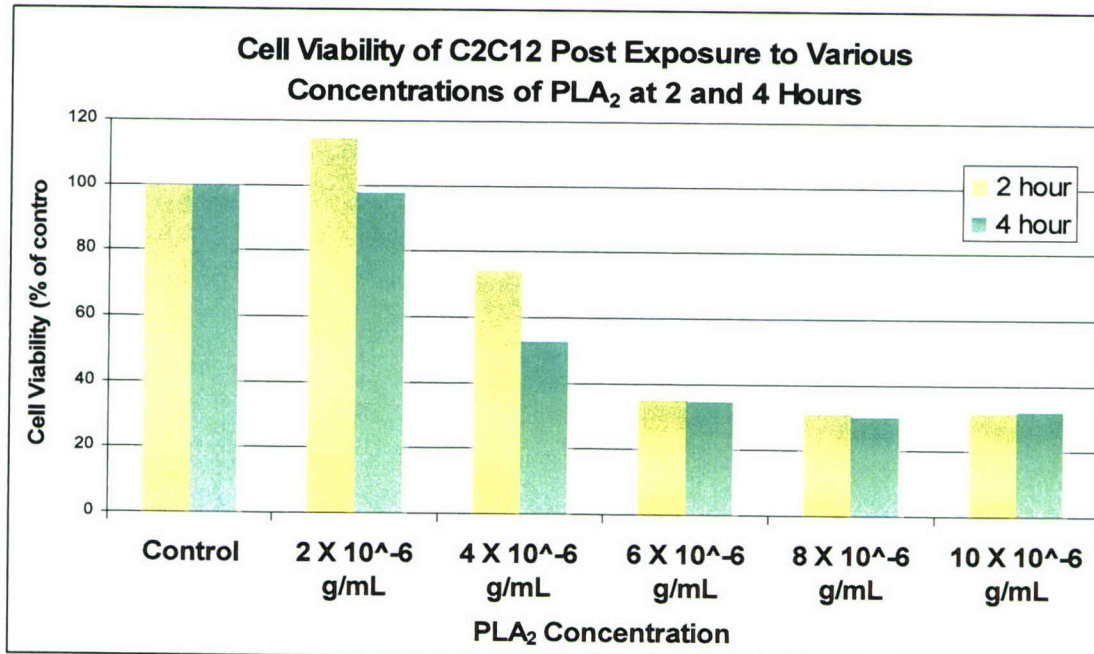


Figure 8. Cell viability of C2C12 cells 2 and 4 hours post exposure to various concentrations of PLA₂.

In the case of MDCK cells, as seen in Figure 9, the LD₅₀ is approximately 30µg/mL for 2 hour exposure periods and approximately 13µg/mL for 4 hour exposure periods.

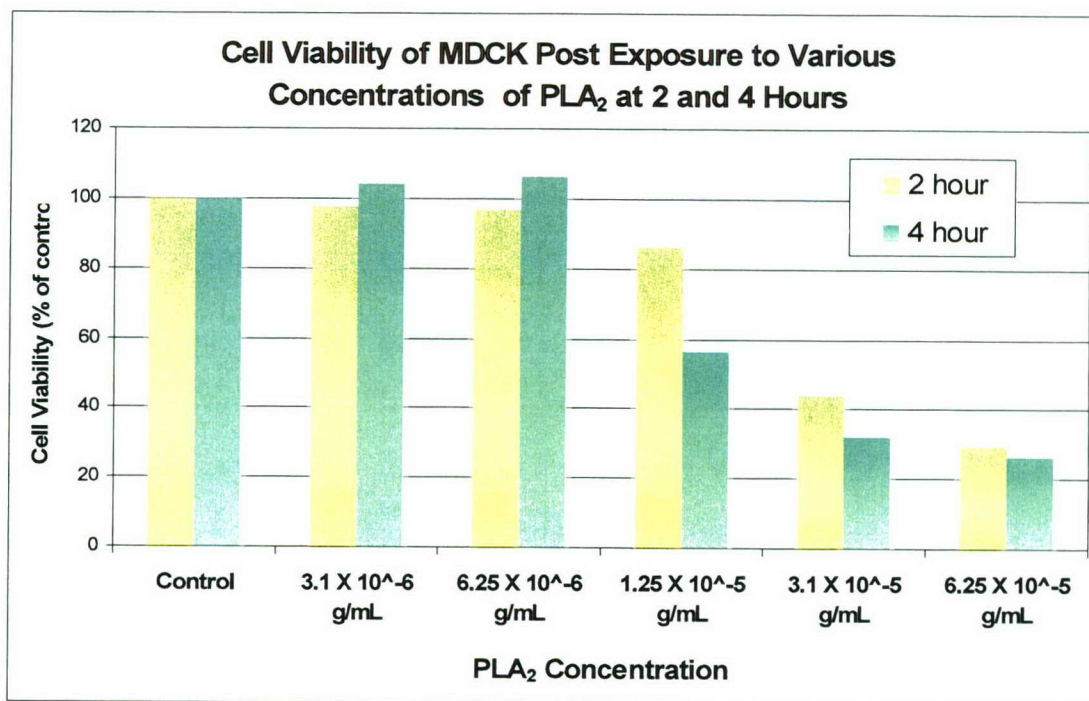


Figure 9. Cell viability of MDCK cells 2 and 4 hours post exposure to various concentrations of PLA₂.

4.3 Inhibitor Evaluation

Inhibitors were tested for their effectiveness against these LD₅₀ values for each cell line. However, in these assays, known LD₅₀ values for PLA₂ did not prove to be toxic to the cells themselves. Higher concentrations of PLA₂, 10µg/mL for C2C12 cells and 50µg/mL for MDCK cells, were also ineffective in killing cells as shown in Figures 10 and 11.

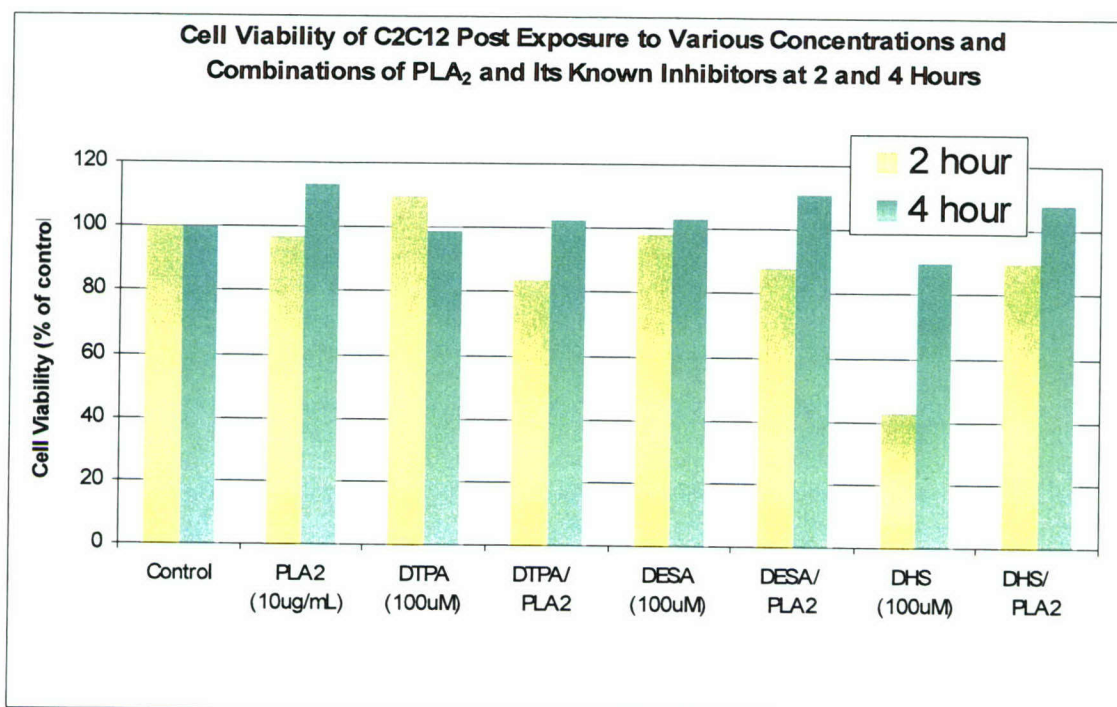


Figure 10. C2C12 viability 2 and 4 hours post exposure to PLA₂ and its known inhibitors.

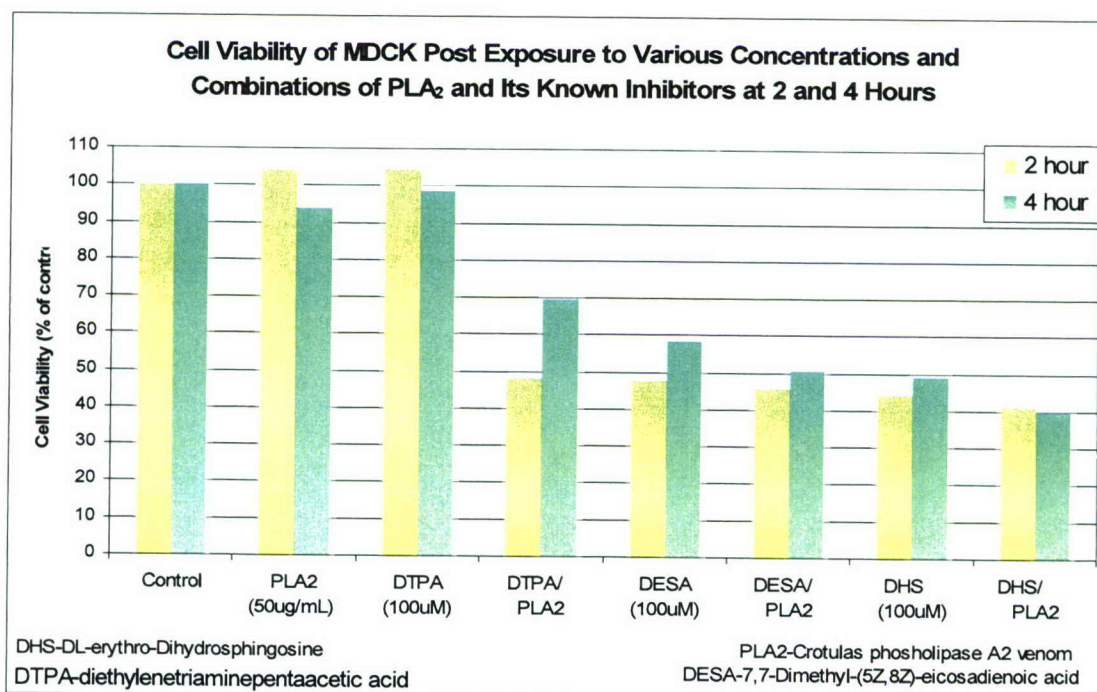


Figure 11. MDCK viability 2 and 4 hours post exposure to PLA₂ and its known inhibitors.

Based on the results presented, the effectiveness of the inhibitor's ability to decrease PLA₂ activity, thereby preventing toxicity to cells, could not be determined.

Chapter 5 MSDS

5.1 PHOSPHOLIPASE A₂ FROM CROTALUS DURISSUS TERRIFICUS

SIGMA-ALDRICH

MATERIAL SAFETY DATA SHEET

Substance Name CAS # SARA 313

PHOSPHOLIPASE A(2) 9001-84-7 No

Synonyms E.C. 3.1.1.4 * Lecithinase A * Phosphatidase * Phosphatidolipase *

Phospholipase A * Phospholipase A(sub 2)

RTECS Number: SZ6114900

Section 3 - Hazards Identification

EMERGENCY OVERVIEW

Poison. May be fatal if enters bloodstream. Do not breathe dust. Do not use if skin is cut or scratched. Wash thoroughly after handling.

HMIS RATING

HEALTH: 4

FLAMMABILITY: 0

REACTIVITY: 0

NFPA RATING

HEALTH: 4

FLAMMABILITY: 0

REACTIVITY: 0

For additional information on toxicity, please refer to Section 11.

Section 4 - First Aid Measures

Section 5 - Fire Fighting Measures

FLASH POINT N/A

AUTOIGNITION TEMP N/A

FLAMMABILITY N/A

EXTINGUISHING MEDIA

Suitable: Water spray. Carbon dioxide, dry chemical powder, or appropriate foam.

FIREFIGHTING

Protective Equipment: Wear self-contained breathing apparatus and protective clothing to prevent contact with skin and eyes. Specific Hazard(s): Emits toxic fumes under fire conditions.

Section 6 - Accidental Release Measures

PROCEDURE TO BE FOLLOWED IN CASE OF LEAK OR SPILL

Evacuate area.

PROCEDURE(S) OF PERSONAL PRECAUTION(S)

Wear self-contained breathing apparatus, rubber boots, and heavy rubber gloves.

METHODS FOR CLEANING UP

Spilled material should be carefully wiped up or moistened with water and removed.

Ventilate area and wash spill site after material pickup is complete.

Section 7 - Handling and Storage

HANDLING

User Exposure: Avoid inhalation. Do not get in eyes, on skin, on clothing. Avoid prolonged or repeated exposure. Do not use if skin is cut or scratched. Wash thoroughly after handling.

STORAGE

Suitable: Keep tightly closed.

Store at -20°C

Section 8 - Exposure Controls / PPE

ENGINEERING CONTROLS

Safety shower and eye bath. Use only in a chemical fume hood.

PERSONAL PROTECTIVE EQUIPMENT

Respiratory: Use respirators and components tested and approved under appropriate government standards such as NIOSH (US) or CEN (EU). Where risk assessment shows air-purifying respirators are appropriate use a full-face particle respirator type N100 (US) or type P3 (EN 143) respirator cartridges as a backup to engineering controls. If the respirator is the sole means of protection, use a full-face supplied air respirator. Hand: Compatible chemical-resistant gloves. Eye: Chemical safety goggles.

GENERAL HYGIENE MEASURES

Wash contaminated clothing before reuse.

Section 9 - Physical/Chemical Properties

Appearance Physical State: Solid Property Value At Temperature or Pressure

pH N/A

BP/BP Range N/A

MP/MP Range N/A

Freezing Point N/A

Vapor Pressure N/A

Vapor Density N/A

Saturated Vapor Conc. N/A

SG/Density N/A

Bulk Density N/A

Odor Threshold N/A
Volatile% N/A
VOC Content N/A
Water Content N/A
Solvent Content N/A
Evaporation Rate N/A
Viscosity N/A
Surface Tension N/A
Partition Coefficient N/A
Decomposition Temp. N/A
Flash Point N/A
Explosion Limits N/A
Flammability N/A
Auto ignition Temp N/A
Refractive Index N/A
Optical Rotation N/A
Miscellaneous Data N/A
Solubility N/A
N/A = not available

Section 10 - Stability and Reactivity

STABILITY Stable: Stable. Materials to Avoid: Strong oxidizing agents.

HAZARDOUS DECOMPOSITION PRODUCTS

Hazardous Decomposition Products: Nature of decomposition products not known.

HAZARDOUS POLYMERIZATION

Hazardous Polymerization: Will not occur

Section 11 - Toxicological Information

ROUTE OF EXPOSURE

Skin Contact: May cause skin irritation.

Skin Absorption: May be harmful if absorbed through the skin.

Eye Contact: May cause eye irritation.

Inhalation: May be harmful if inhaled. Material may be irritating to mucous membranes and upper respiratory tract.

Ingestion: May be harmful if swallowed.

SENSITIZATION

Sensitization: Prolonged or repeated exposure may cause allergic reactions in certain sensitive individuals.

SIGNS AND SYMPTOMS OF EXPOSURE

To the best of our knowledge, the chemical, physical, and toxicological properties have not been thoroughly investigated.

CONDITIONS AGGRAVATED BY EXPOSURE

Basic phospholipases in general appear to be toxic, having lethal dose of about 500 ug/kg/mouse or even considerably less as exemplified by presynaptic neurotoxins and

myonecrotins. Many of the acidic phospholipases A, however, are much less toxic and not lethal even at 2000 ug/kg, although there are exceptions. An acidic phospholipase a (isoelectric point 5.1) from naja nigricollis has a lethal dose of 800 ug/kg in the mouse.* The toxicological properties have not been thoroughly investigated. May be fatal if enters bloodstream.

TOXICITY DATA

Intraperitoneal Mouse 5 MG/KG LD50
Subcutaneous Mouse 1200 UG/KG LD50
Intravenous Mouse 7500 UG/KG LD50

Section 12 - Ecological Information

No data available.

Section 13 - Disposal Considerations

APPROPRIATE METHOD OF DISPOSAL OF SUBSTANCE OR PREPARATION
Contact a licensed professional waste disposal service to dispose of this material.
Observe all federal, state, and local environmental regulations.

Section 14 - Transport Information DOT

Proper Shipping Name: None
Non-Hazardous for Transport: This substance is considered to be non-hazardous for transport.
IATA
Non-Hazardous for Air Transport: Non-hazardous for air transport.

Section 15 - Regulatory Information

US CLASSIFICATION AND LABEL TEXT
US Statements: Poison. May be fatal if enters bloodstream. Do not breathe dust. Do not use if skin is cut or scratched. Wash thoroughly after handling.
UNITED STATES REGULATORY INFORMATION
SARA LISTED: No
CANADA REGULATORY INFORMATION
WHMIS Classification: This product has been classified in accordance with the hazard criteria of the CPR, and the MSDS contains all the information required by the CPR.
DSL: No
NDSL: No

Section 16 - Other Information

DISCLAIMER

For R&D use only. Not for drug, household or other uses.

5.2 DIETHYLENETRIAMINEPENTAACETIC ACID

SIGMA-ALDRICH

MATERIAL SAFETY DATA SHEET

Section 1 - Product and Company Information

Product Name DIETHYLENETRIAMINEPENTAACETIC ACID FREE&

Product Number D6518

Brand SIAL

Company Sigma-Aldrich

Address 3050 Spruce Street

SAINT LOUIS MO 63103 US

Technical Phone: 800-325-5832

Fax: 800-325-5052

Emergency Phone: 314-776-6555

Section 2 - Composition/Information on Ingredient

Substance Name CAS # SARA 313

DIETHYLENETRIAMINE PENTAACETIC ACID 67-43-6 Yes

Ingredient Name CAS # Percent SARA 313

NITRILOTRIACETIC ACID 139-13-9 0.2 Yes

Formula C₁₄H₂₃N₃O₁₀

Synonyms Acetic acid,

((carboxymethylimino)bis(ethylenenitrilo))tetra-*

((((Carboxymethyl)imino)bis(ethylenenitrilo))tetraa

cetic acid * CHEL 330 * CHEL 330 acid * Chel DTPA

* Dabeersen 503 * Detapac * Detarex *

Diethylenetriaminepentaacetic acid *

1,1,4,7,7-Diethylenetriaminepentaacetic acid *

(Diethylenetrinitrilo)pentaacetic acid * DTPA * Hamp-EX Acid * Monaquest CAI *

Pentetic acid * Titriplex V * 3,6,9-Triazaundecanedioic acid, 3,6,9-tris(carboxymethyl)-

RTECS Number: MB8205000

Section 3 - Hazards Identification

EMERGENCY OVERVIEW

Irritant. Dangerous for the environment. Irritating to eyes. Toxic to aquatic organisms, may cause long-term adverse effects in the aquatic environment.

Calif. Prop. 65 carcinogen.

HMIS RATING

HEALTH: 2*

FLAMMABILITY: 0

REACTIVITY: 0

NFPA RATING

HEALTH: 2

FLAMMABILITY: 0

REACTIVITY: 0

*additional chronic hazards present.

For additional information on toxicity, please refer to Section 11.

Section 4 - First Aid Measures

ORAL EXPOSURE

If swallowed, wash out mouth with water provided person is conscious. Call a physician.

INHALATION EXPOSURE

If inhaled, remove to fresh air. If not breathing give artificial respiration. If breathing is difficult, give oxygen.

DERMAL EXPOSURE

In case of contact, immediately wash skin with soap and copious amounts of water.

EYE EXPOSURE

In case of contact, immediately flush eyes with copious amounts of water for at least 15 minutes.

Section 5 - Fire Fighting Measures

FLASH POINT

392 °F 200 °C Method: closed cup

AUTOIGNITION TEMP N/A

FLAMMABILITY N/A

EXTINGUISHING MEDIA

Suitable: Water spray. Carbon dioxide, dry chemical powder, or appropriate foam.

FIREFIGHTING

Protective Equipment: Wear self-contained breathing apparatus and protective clothing to prevent contact with skin and eyes.

Specific Hazard(s): Emits toxic fumes under fire conditions.

Section 6 - Accidental Release Measures

PROCEDURE(S) OF PERSONAL PRECAUTION(S)

Wear respirator, chemical safety goggles, rubber boots, and heavy rubber gloves.

METHODS FOR CLEANING UP

Sweep up, place in a bag and hold for waste disposal. Avoid raising dust. Ventilate area and wash spill site after material pickup is complete.

Section 7 - Handling and Storage

HANDLING

User Exposure: Do not breathe dust. Avoid contact with eyes, skin, and clothing. Avoid prolonged or repeated exposure.

STORAGE

Suitable: Keep tightly closed.

Section 8 - Exposure Controls / PPE

ENGINEERING CONTROLS

Safety shower and eye bath. Mechanical exhaust required.

PERSONAL PROTECTIVE EQUIPMENT

Respiratory: Use respirators and components tested and approved under appropriate government standards such as NIOSH (US) or CEN (EU). Where risk assessment shows air-purifying respirators are appropriate use a dust mask type N95 (US) or type P1 (EN 143) respirator.

Hand: Compatible chemical-resistant gloves.

Eye: Chemical safety goggles.

GENERAL HYGIENE MEASURES

Wash thoroughly after handling.

Section 9 - Physical/Chemical Properties

Appearance Physical State: Solid

Property Value At Temperature or Pressure

Molecular Weight 393.35 AMU

pH 2.5 23 °C Concentration: 10g/l

BP/BP Range N/A

MP/MP Range 219.0 - 220.0 °C

Freezing Point N/A

Vapor Pressure N/A

Vapor Density N/A

Saturated Vapor Conc. N/A

SG/Density N/A

Bulk Density N/A

Odor Threshold N/A

Volatile% N/A

VOC Content N/A

Water Content N/A

Solvent Content N/A

Evaporation Rate N/A

Viscosity N/A

Surface Tension N/A

Partition Coefficient N/A

Decomposition Temp. N/A

Flash Point 392 °F 200 °C Method: closed cup

Explosion Limits N/A

Flammability N/A

Auto ignition Temp N/A

Refractive Index N/A

Optical Rotation N/A
Miscellaneous Data N/A
Solubility in Water: 1 g/l, 20°C
Solvent: clear, colorless 0.1M in NaOH 1M,
20°C
N/A = not available

Section 10 - Stability and Reactivity

STABILITY

Stable: Stable.

Materials to Avoid: Strong oxidizing agents.

HAZARDOUS DECOMPOSITION PRODUCTS

Hazardous Decomposition Products: Carbon monoxide, Carbon dioxide,
Nitrogen oxides.

HAZARDOUS POLYMERIZATION

Hazardous Polymerization: Will not occur

Section 11 - Toxicological Information

ROUTE OF EXPOSURE

Skin Contact: May cause skin irritation.

Skin Absorption: May be harmful if absorbed through the skin.

Eye Contact: Causes eye irritation.

Inhalation: Material may be irritating to mucous membranes and upper respiratory tract.
May be harmful if inhaled.

Ingestion: May be harmful if swallowed.

SIGNS AND SYMPTOMS OF EXPOSURE

To the best of our knowledge, the chemical, physical, and toxicological properties have
not been thoroughly investigated.

TOXICITY DATA

Oral Rat > 2,000 mg/kg LD50

Intraperitoneal Rat 587 MG/KG LD50

Remarks: Behavioral: Convulsions or effect on seizure threshold.

Behavioral: Aggression. Lungs, Thorax, or Respiration: Chronic
pulmonary edema.

Intraperitoneal Mouse 543 MG/KG LD50

IRRITATION DATA

Skin Rabbit

Remarks: No irritation effect

Eyes Rabbit

Remarks: Moderate irritation effect

Section 12 - Ecological Information

ACCUMULATION

Bioaccumulation Potential: Indication of bioaccumulation.

ACUTE ECOTOXICITY TESTS

Test Type: flow-through bioassay

Species: *Leuciscus idus*

Time: 96 h

Value: > 100 mg/l

SIALL - D6518 www.sigma-aldrich.com Page 4

Test Type: Growth inhibitor on algae.

Time: 72 h

Value: 1.0 - 10.0 mg/l

Test Type: EC50 *Daphnia*

Species: *Daphnia*

Time: 48 h

Value: 245 mg/l

ELIMINATION

Elimination: 20.0 - 60.0 %

Section 13 - Disposal Considerations

APPROPRIATE METHOD OF DISPOSAL OF SUBSTANCE OR PREPARATION

Contact a licensed professional waste disposal service to dispose of this material.

Dissolve or mix the material with a combustible solvent and burn in a chemical incinerator equipped with an afterburner and scrubber. Observe all federal, state, and local environmental regulations.

Section 14 - Transport Information

DOT

Proper Shipping Name: Environmentally hazardous substances, solid, n.o.s.

UN#: 3077

Class: 9

Packing Group: Packing Group III

Hazard Label: Class 9

PIH: Not PIH

IATA

Proper Shipping Name: Environmentally hazardous substance, solid, n.o.s

IATA UN Number: 3077

Hazard Class: 9

Packing Group: III

Section 15 - Regulatory Information

EU ADDITIONAL CLASSIFICATION

Symbol of Danger: Xi-N

Indication of Danger: Irritant. Dangerous for the environment.

R: 36-51/53

Risk Statements: Irritating to eyes. Toxic to aquatic organisms, may cause long-term adverse effects in the aquatic environment.

S: 26-36-61

Safety Statements: In case of contact with eyes, rinse immediately with plenty of water and seek medical advice. Wear suitable protective clothing. Avoid release to the environment. Refer to special instructions/safety data sheets.

US CLASSIFICATION AND LABEL TEXT

Indication of Danger: Irritant. Dangerous for the environment.

Risk Statements: Irritating to eyes. Toxic to aquatic organisms, may cause long-term adverse effects in the aquatic environment.

Safety Statements: In case of contact with eyes, rinse immediately with plenty of water and seek medical advice. Wear suitable protective clothing. Avoid release to the environment.

Refer to special instructions/safety data sheets.

US Statements: Calif. Prop. 65 carcinogen.

UNITED STATES REGULATORY INFORMATION

SARA LISTED: Yes

NOTES: This product is or contains a component that is subject to SARA313 reporting requirements.

TSCA INVENTORY ITEM: Yes

UNITED STATES - STATE REGULATORY INFORMATION

CALIFORNIA PROP - 65

California Prop - 65: This product is or contains chemical(s) known to the state of California to cause cancer.

CANADA REGULATORY INFORMATION

WHMIS Classification: This product has been classified in accordance with the hazard criteria of the CPR, and the MSDS contains all the information required by the CPR.

DSL: Yes

NDL: No

Section 16 - Other Information

DISCLAIMER

For R&D use only. Not for drug, household or other uses.

5.3 DL-ERYTHRO-DIHYDROSPHINGOSINE

SIGMA-ALDRICH

MATERIAL SAFETY DATA SHEET

Section 1 - Product and Company Information

Product Name DL-ERYTHRO-DIHYDROSPHINGOSINE

Product Number D6908
Brand SIGMA
Company Sigma-Aldrich
Address 3050 Spruce Street
SAINT LOUIS MO 63103 US
Technical Phone: 800-325-5832
Fax: 800-325-5052
Emergency Phone: 314-776-6555

Section 2 - Composition/Information on Ingredient

Substance Name CAS # SARA 313
DL-ERYTHRO-DIHYDROSPHINGOSINE 3102-56-5 No
Formula $C_{18}H_{39}NO_2$

Section 3 - Hazards Identification

HMIS RATING

HEALTH: 0

FLAMMABILITY: 0

REACTIVITY: 0

NFPA RATING

HEALTH: 0

FLAMMABILITY: 0

REACTIVITY: 0

For additional information on toxicity, please refer to Section 11.

Section 4 - First Aid Measures

Section 5 - Fire Fighting Measures

FLASH POINT N/A

AUTOIGNITION TEMP N/A

FLAMMABILITY N/A

Section 6 - Accidental Release Measures

Section 7 - Handling and Storage

STORAGE Store at -20°C

Section 8 - Exposure Controls / PPE

Section 9 - Physical/Chemical Properties

Appearance Color: White

Form: Powder
Property Value At Temperature or Pressure
Molecular Weight 301.5 AMU
pH N/A
BP/BP Range N/A
MP/MP Range N/A
Freezing Point N/A
Vapor Pressure N/A
Vapor Density N/A
Saturated Vapor Conc. N/A
SG/Density N/A
Bulk Density N/A
Odor Threshold N/A
Volatile% N/A
VOC Content N/A
Water Content N/A
Solvent Content N/A
Evaporation Rate N/A
Viscosity N/A
Surface Tension N/A
Partition Coefficient N/A
Decomposition Temp. N/A
Flash Point N/A
Explosion Limits N/A
Flammability N/A
Auto ignition Temp N/A
Refractive Index N/A
Optical Rotation N/A
Miscellaneous Data N/A
Solubility Solvent: clear, colorless 20 mg/ml CHCl₃
N/A = not available

Section 10 - Stability and Reactivity

Section 11 - Toxicological Information

Section 12 - Ecological Information

No data available.

Section 13 - Disposal Considerations

Section 14 - Transport Information

DOT
Proper Shipping Name: None

Non-Hazardous for Transport: This substance is considered to be non-hazardous for transport.

IATA

Non-Hazardous for Air Transport: Non-hazardous for air transport.

Section 15 - Regulatory Information

UNITED STATES REGULATORY INFORMATION

SARA LISTED: No

CANADA REGULATORY INFORMATION

WHMIS Classification: This product has been classified in accordance with the hazard criteria of the CPR, and the MSDS contains all the information required by the CPR.

DSL: No

NDSL: No

Section 16 - Other Information

DISCLAIMER

For R&D use only. Not for drug, household or other uses.

5.4 7,7-DIMETHYL-(5Z,8Z)-EICOSADIENOIC ACID

SIGMA-ALDRICH

MATERIAL SAFETY DATA SHEET

Section 1 - Product and Company Information

Product Name 7,7-DIMETHYL-(5Z,8Z)-EICOSADIENOIC ACID

Product Number D8008

Brand SIGMA

Company Sigma-Aldrich

Address 3050 Spruce Street

SAINT LOUIS MO 63103 US

Technical Phone: 800-325-5832

Fax: 800-325-5052

Emergency Phone: 314-776-6555

Section 2 - Composition/Information on Ingredient

Substance Name CAS # SARA 313

7,7-DIMETHYL-(5Z,8Z)-EICOSADIENOIC 89560-01-0 Yes
ACID

Ingredient Name CAS # Percent SARA 313

The hazards identified with this None product are those associated with the residual solvents that are used in the manufacturing process.

METHANOL 67-56-1 <= 1 Yes

Formula C₂₂H₄₀O₂

Section 3 - Hazards Identification

EMERGENCY OVERVIEW

Toxic. Toxic by inhalation, in contact with skin and if swallowed. Toxic: danger of very serious irreversible effects through inhalation, in contact with skin and if swallowed. Irritating to eyes and skin.

Target organ(s): Eyes. Kidneys.

HMIS RATING

HEALTH: 2*

FLAMMABILITY: 0

REACTIVITY: 1

NFPA RATING

HEALTH: 2

FLAMMABILITY: 0

REACTIVITY: 1

*additional chronic hazards present.

For additional information on toxicity, please refer to Section 11.

Section 4 - First Aid Measures

ORAL EXPOSURE

If swallowed, wash out mouth with water provided person is conscious. Call a physician immediately.

INHALATION EXPOSURE

If inhaled, remove to fresh air. If not breathing give artificial respiration. If breathing is difficult, give oxygen.

DERMAL EXPOSURE

In case of skin contact, flush with copious amounts of water for at least 15 minutes. Remove contaminated clothing and shoes. Call a physician.

EYE EXPOSURE

In case of contact with eyes, flush with copious amounts of water for at least 15 minutes. Assure adequate flushing by separating the eyelids with fingers. Call a physician.

Section 5 - Fire Fighting Measures

FLASH POINT N/A

AUTOIGNITION TEMP N/A

FLAMMABILITY N/A

EXTINGUISHING MEDIA

Suitable: Water spray. Carbon dioxide, dry chemical powder, or appropriate foam.

FIREFIGHTING

Protective Equipment: Wear self-contained breathing apparatus and protective clothing to prevent contact with skin and eyes.

Section 6 - Accidental Release Measures

PROCEDURE(S) OF PERSONAL PRECAUTION(S)

Wear respirator, chemical safety goggles, rubber boots, and heavy rubber gloves.

METHODS FOR CLEANING UP

Absorb on sand or vermiculite and place in closed containers for disposal. Ventilate area and wash spill site after material pickup is complete.

Section 7 - Handling and Storage

HANDLING

User Exposure: Do not breathe vapor. Avoid contact with eyes, skin, and clothing. Avoid prolonged or repeated exposure.

STORAGE

Suitable: Keep container closed.

Store at -20°C

SPECIAL REQUIREMENTS

Light sensitive.

Section 8 - Exposure Controls / PPE

ENGINEERING CONTROLS

Safety shower and eye bath. Mechanical exhaust required.

PERSONAL PROTECTIVE EQUIPMENT

Respiratory: Use respirators and components tested and approved under appropriate government standards such as NIOSH (US) or CEN (EU). Use supplied-air or SCBA respirators. Europe permits the use of type AXBEK full-face cartridge respirators (EN 14387).

Other: Wear appropriate government approved respirator, chemical-resistant gloves, safety goggles, other protective clothing.

GENERAL HYGIENE MEASURES

Wash thoroughly after handling.

EXPOSURE LIMITS

Country Source Type Value

Poland NDS 100 MG/M3

Poland NDSCh 300 MG/M3

Poland NDSP –

Section 9 - Physical/Chemical Properties

Appearance Physical State: Liquid

Property Value At Temperature or Pressure

Molecular Weight 336.6 AMU

pH N/A

BP/BP Range N/A

MP/MP Range N/A

Freezing Point N/A

Vapor Pressure N/A

Vapor Density N/A

Saturated Vapor Conc. N/A

SG/Density N/A

Bulk Density N/A

Odor Threshold N/A

Volatile% N/A

VOC Content N/A
Water Content N/A
Solvent Content N/A
Evaporation Rate N/A
Viscosity N/A
Surface Tension N/A
Partition Coefficient N/A
Decomposition Temp. N/A
Flash Point N/A
Explosion Limits N/A
Flammability N/A
Auto ignition Temp N/A
Refractive Index N/A
Optical Rotation N/A
Miscellaneous Data N/A
Solubility N/A
N/A = not available

Section 10 - Stability and Reactivity

STABILITY

Stable: Stable.

Conditions to Avoid: Exposure to light may affect product quality.

Materials to Avoid: Strong oxidizing agents.

HAZARDOUS DECOMPOSITION PRODUCTS

Hazardous Decomposition Products: Carbon monoxide, Carbon dioxide.

HAZARDOUS POLYMERIZATION

Hazardous Polymerization: Will not occur

Section 11 - Toxicological Information

ROUTE OF EXPOSURE

Skin Contact: May cause skin irritation.

Skin Absorption: May be harmful if absorbed through the skin.

Eye Contact: May cause eye irritation.

Inhalation: Material may be irritating to mucous membranes and upper respiratory tract.

May be harmful if inhaled.

Ingestion: May be harmful if swallowed.

TARGET ORGAN(S) OR SYSTEM(S)

Eyes. Kidneys. Liver. Heart.

SIGNS AND SYMPTOMS OF EXPOSURE

To the best of our knowledge, the chemical, physical, and toxicological properties have not been thoroughly investigated. May cause convulsions. Gastrointestinal disturbances.

Section 12 - Ecological Information

No data available.

Section 13 - Disposal Considerations

APPROPRIATE METHOD OF DISPOSAL OF SUBSTANCE OR PREPARATION

Contact a licensed professional waste disposal service to dispose of this material. Dissolve or mix the material with a combustible solvent and burn in a chemical incinerator equipped with an afterburner and scrubber. Observe all federal, state, and local environmental regulations.

Section 14 - Transport Information

DOT

Proper Shipping Name: None

Non-Hazardous for Transport: This substance is considered to be non-hazardous for transport.

IATA

Non-Hazardous for Air Transport: Non-hazardous for air transport.

Section 15 - Regulatory Information

US CLASSIFICATION AND LABEL TEXT

Indication of Danger: Toxic.

Risk Statements: Toxic by inhalation, in contact with skin and if swallowed. Toxic: danger of very serious irreversible effects through inhalation, in contact with skin and if swallowed. Irritating to eyes and skin.

Safety Statements: Keep container tightly closed. Avoid contact with skin. Wear suitable protective clothing and gloves. In case of accident or if you feel unwell, seek medical advice immediately (show the label where possible).

US Statements: Target organ(s): Eyes. Kidneys.

UNITED STATES REGULATORY INFORMATION

SARA LISTED: Yes

NOTES: This product is or contains a component that is subject to SARA313 reporting requirements.

CANADA REGULATORY INFORMATION

WHMIS Classification: This product has been classified in accordance with the hazard criteria of the CPR, and the MSDS contains all the information required by the CPR.

DSL: No

NDSL: No

Section 16 - Other Information

DISCLAIMER

For R&D use only. Not for drug, household or other uses.

All information, recommendations and suggestions herein concerning this product are based upon data believed to be reliable. However it is the user's responsibility to determine the safety, toxicity and suitability for his/her own use of this product. Since the actual use of others is beyond our control, we make no guarantee expressed or implied as to the effects of such use, the results to be obtained, or the safety and toxicity of the product. This information is not to be construed as absolutely complete, since additional information may be necessary of desirable when exceptional conditions or circumstances exist or because of applicable laws or government regulations.

Disclaimer (provided with this information by the compiling agencies): This information is formulated for use by elements of the Department of Defense. The United States of America in no manner whatsoever expressly or implied warrants, states, or intends said information to have any application, use or viability by or to any person or persons outside the Department of Defense nor any person or persons contracting with any instrumentality of the United States of America and disclaims all liability for such use. Any person utilizing this instruction who is not a military or civilian employee of the United States of America should seek competent professional advice to verify and assume responsibility for the suitability of this information to their particular situation regardless of similarity to a corresponding Department of Defense or other government situation.

Chapter 6 References

Reference List

1. Andriao-Escarso, Silvia, Soares, Andreimar, Rodrigues, Veidiana, Angulo, Yamileth, Diaz, Cecilia, Lomonte, Bruno, Gutierrez, Jose and Giglio, Jose. *Biochimie* 82, 755-63, (2000).
2. Soares, Andreimar, Marcussi, Silvana, Stabeli, Rodrigo, Franca, Suzelei, Giglio, Jose, Ward, Richard, and Arantes, Eliane. *Biochem Biophys Res Commun* 302, 193-200 (2003).
3. Floege, Jurgen, Ostendorf, Tammo, Janssen, Ulf, Burg, Michael, Radeke, Heinfried, Vargeese, Chandra, Gill, Stanley, Green, Louis, Janjic, Nebojsa. *Am.J.Pathol.* 154, 169-179 (1999).
4. Paddle, Brian. *J Appl Toxicol.* 23, 139-170 (2003).
5. Brocker, Fred, Hanson, Kevin, and Adams, Ann. Campaign Analysis Report: Venomous Snakes and Scorpions in Iraq, and Their Antivenin Sources.
<<http://www.brooks.af.mil/web/af/courses/amp/cluebag/Venomous%20Snakes%20&%20Scorpions%20in%20Iraq%20&%20Antivenin%20Sources.pdf>> U-145,017-02. 8-10-2002. Joint Chiefs of Staff . 8-10-2002.
Ref Type: Electronic Citation
6. Radmanesh, Mohammed. *Inter J Dermatol* 37, 500-507 (1998).
7. Radmanesh, Mohammed. *J Trop.Med Hyg.* 93, 327-332 (1990).
8. Chadha, J.S. and Leviav, A. *JAMA* 241, 1038 (1979).
9. Bradbury, Andrew and Cattaneo, Antonio. *Trends Neurosci.* 18, 243-249 (1995).
10. Clackson, Tim and Wells, James. *Trends Biotechnol.* 12, 173-184 (1994).
11. Tuerk, Craig and Gold, Larry *Science* 249, 505-510 (1990).
12. Ellington, Andrew and Szostak, Jack. *Nature* 346, 818-822 (1990).
13. Ciesiolka, Jerzy and Yarus, Michael. *RNA.* 2, 785-793 (1996).
14. Bless, Nicolas, Smith, Drew, Charlton, Josephine, Czermak, Boris, Schmal, Hagen, Friedl, Hans, and Ward, Peter, *Curr.Biol.* 7, 877-880 (1997).

15. Lebruska, Lori and Maher III, L. James, *Biochemistry* 38, 3168-3174 (1999).
16. Hesselberth, Jay, Miller, Darcie, Robertus, Jon, and Ellington, Andrew. *J.Biol.Chem.* 275, 4937-4942 (2000).
17. Kiel, Johnathan L., Parker, Jill E., Holwitt, Eric A., and Vivekananda, Jeeva. DNA capture elements for rapid detection and identification of biological agents. Gardner, Patrick J. 5416(1), 105-110. 8-13-2004. Orlando, FL, USA, SPIE.
Ref Type: Conference Proceeding
18. Li, Jianwei, Fang, Xiaohong, and Tan, Weihong. *Biochem Biophys Res Commun* 292, 31-40 (2002).
19. Kensch, Oliver, Connolly, Bernard, Steinhoff, Heinz-Jurgen, McGregor, Alistair, Goody, Roger, and Restle, Tobias *J.Biol Chem.* 275, 18271-8 (2000).
20. Zhang, Xing-Mei, Li, Quian., Shi, Yu-Sheng, Wu, Jun-Hua, Shao, Ning-Sheng, Liu, Gang, Sun, Man-Ji. *Brain Res.* 989, 147-151 (2003).
21. Vinores, S.A. *Curr.Opin.Mol.Ther.* 5, 673-679 (2003).
22. Osborne, Scott and Ellington, Andrew. *Chem Rev.* 97, 349-370 (1997).
23. Vivekananda, Jeeva and Kiel, Johnathan L. Methods and compositions for aptamers against anthrax. Conceptual MindWorks, Inc. 978753(6569630). 2005. Texas/United States of America.
Ref Type: Patent
24. Smith, Theodore. and Figge, Helen. *Am J Health Syst Pharm.* 48, 2190-2196 (1991).
25. Ruckman, Judy, Green, Louis, Beeson, Jim, Waugh, Sheela, Gillette, Wendy, Henninger, Dwight, Claesson-Welsh, Lena, and Janjic, Neboja. *J Biol.Chem* 273, 20556-20567 (1998).
26. Tucker, Chistopher, Chen, Long-Shiuh, Judkins, Mark, Farmer, James, Gill, Stanley, and Drolet, Daniel. *J Chromatogr.B Biomed.Sci Appl* 732, 203-212 (1999).
27. Gerlier, Denise and Thomasset, Nicole. *J Immunol Methods* 94, 57-63 (1986)